

Material Compatibility List

for Exact Imaging ExactVu™ Transducers

Material Compatibility List (MCL): Listed Chemicals do not degrade transducer performance or safety over 2,500 reprocessing cycles (500 cycles for abdominal high level disinfectants).

NOTE: Microbiological efficacy of the chemicals on the MCL has not been evaluated by Exact Imaging. Please refer to chemical manufacturers' claims.

Type of Chemical	Chemical Name	Manufacturer	Material Compatibility			
			ExactVu™ Transducers			
			EV9C	EV9C	EV29L	EV5C
			SN starting 6087-	SN starting 6750-	SN starting 6089-	SN starting 6088-
Cleaners	Cidezyme	Advanced Sterilization Products	✓	✓	✓	✓
	Enzol	Advanced Sterilization Products	✓	✓	✓	✓
	Aniosyme X3	Laboratories Anios	NC	✓	✓	✓
	Aniosyme DD1	Laboratories Anios	NC	✓	✓	✓
	Revital-Ox Cleaner	Steris	NC	✓	✓	✓
	Transeptic Spray	Parker Labs	NT	NT	NT	✓

High Level Disinfectants	Aldahol 1.8	Olympus America	NT	NT	✓	NT
	Anioxyde 1000	Laboratories Anios	NC	✓	✓	✓
	Cidex Activated Dialdehyde	Advanced Sterilization Products	✓	✓	✓	✓
	Cidex OPA	Advanced Sterilization Products	✓	✓	✓	✓
	Metricide Activated Dialdehyde	Metrex Research Inc	✓	✓	✓	✓
	Revital-Ox Resert	Steris	NC	✓	✓	✓
	RelyOn PeraSafe	Dupont	NT	NT	✓	NT
	UltrOx	Civco	NC	✓	✓	✓
Wipes/Sprays	Acryl-des	Schulke	✓	✓	✓	✓
	Cleanisept Wipes	Dr. Schumacher GmbH	✓	✓	✓	✓
	Cleanisept Wipes Forte	Dr. Schumacher GmbH	✓	✓	✓	✓
	Incidin Foam	Ecolab	NT	✓	✓	✓
	Tristel Duo	Tristel	✓	✓	✓	✓
	Tristel Trio	Tristel	✓	✓	✓	✓
	Tristel ULT	Tristel	✓	✓	✓	✓
	Surfa'safe	Laboratories Anios	NT	✓	✓	✓
	Wip'Anios	Laboratories Anios	NT	✓	✓	✓

	Clinell Universal Wipes	GAMA Healthcare	NT	✓	✓	✓
	Mikrobac® Tissues	BODE Chemie GmbH	NT	✓	✓	✓
	mikrozid® universal wipes	schülke	NT	✓	✓	✓
	mikrozid® sensitive wipes premium	schülke	NT	✓	✓	✓

Legend:

NC = Not Compatible

NT = Not Tested

Soaking Instructions for ExactVu Transducers:

The images below show the soak level for the transducers.

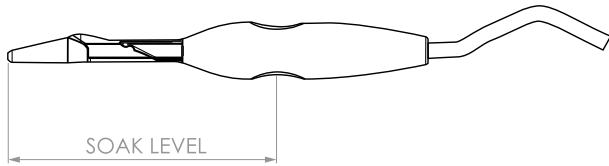


Figure 1: EV29L Soak Level

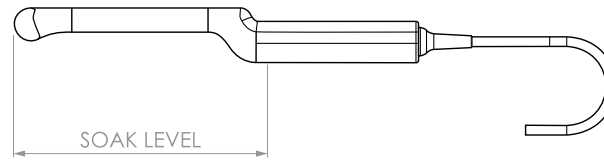


Figure 2: EV9C Soak Level

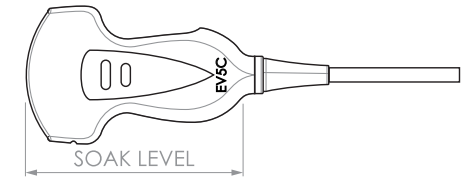


Figure 3: EV5C Soak Level

CAUTION
EN-C25



Do not immerse the transducer beyond its soak level.

Automatic Reprocessor Compatibility

NOTE: For automatic reprocessors please consult Exact Imaging prior to use. Soak levels for automatic reprocessors vary with device, it is important to review with Exact Imaging prior to use to ensure that soak levels higher than pictured in the diagram below do not damage your transducers.

Automatic Reprocessor Name	Manufacturer	Material Compatibility			
		ExactVu™ Transducers			
		EV9C	EV9C	EV29L	EV5C
		SN starting 6087-	SN starting 6750-	SN starting 6089-	SN starting 6088-

ASTRA VR	CIVCO	NT	NT	✓	NT
Advantage Plus	Medivators	NT	NT	✓	NT
DSD Edge	Medivators	NT	NT	✓	NT
DSD-201	Medivators	NT	NT	✓	NT

Legend:

NC = Not Compatible

NT = Not Tested



Refer to IFU (Care, Cleaning, and Use Guide for EV29L/EV9C/EV5C) for further instructions

Over time, cleaning and disinfection may cause discoloration.
Discoloration does not affect the performance of the transducer.

For any further questions about *Cleaning and Disinfection*
with the Exact Imaging transducers,
for North American customers, please contact Apps@edap-usa.com
For Rest of the World customers, please contact IA@edap-tms.com